




Superior optics, high magnification and wide zoom range designed in industrial grade aluminium housing for Full High Definition image quality, where and how you need it! Available in two industrial grade form factors, F30 with integrated stand and F30s modular.



- Superior 30:1 zoom optics with auto-focus
- On board easy-to-access lens- and camera controls
- Large practical working distance
- Fully flexible and configurable with different stands and mounts
- Full HD 1080p, 60 fps for lag-free magnified working
- Options include ESD-protected housing and a built-in laser pointer

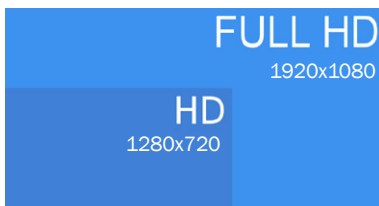

"Efficient, Ergonomic, endlessly versatile!"

Designed for ergonomic use


Inspectis F30s provides relief to eyes, neck and shoulders through the unique ergonomic design, built-in laser pointer and 230 mm free working distance. By allowing operators to sit comfortably in a good working position, their inspection tasks can be carried out more efficiently.




Thanks to its fast optics and large image sensor, F30s provides an excellent image of inspecting object, especially solder joints on electronic PCBs, highly reflective metallic parts and assemblies with minimum glare and reflections.



The system's standard optics with 1.9 - 56x magnification (24" monitor) are comparable to common inspection microscopes. If further magnification is needed, a wide range of +1, +2, +3, +5 and +10 diopter Macro Lenses can be added to magnify up to 170x.

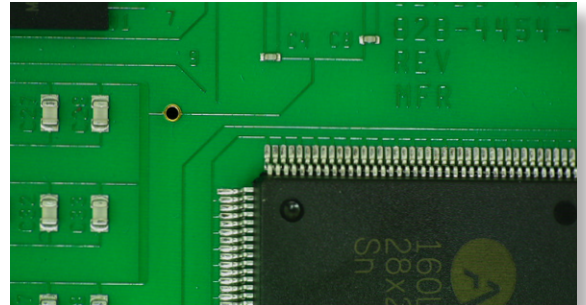

Electronics PCB, x6 Zoom, 11.5x Mag.

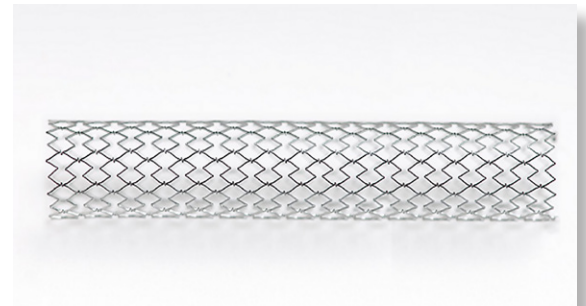

Medtech, x30 Zoom with +5 Macro, 110x Mag.



Micro mechanism, x10 Zoom, 19x Mag.



Forensic, x25 Zoom, 48x Mag.

Connect F30s digital microscope to a Full HD display or your PC via Inspectis' plug-&-play USB3.0 converter and take advantage of the astonishing image resolution, high contrast, true colours at 1:30x zoom range in your optical inspection work.

"Optimized for inspection and repair of electronic PCBs"

Making inspection reliable

F30 offers true ease of use with minimal interaction needed with controls and settings. All basic parameters such as zoom, brightness and colour levels can be controlled by the user through on-board buttons on top of the device.

When required, the optional Advanced Remote Control allows the excellence of the imaging system to be exploited further and additional useful inspection aids to be selected.

To further simplify operations, selected controller functions can be controlled by a foot pedal.



- ☐ 1:30x Zoom Optics
- ☐ AutoFocus
- ☐ On board Controls
- ☐ Aluminium Housing
- ☐ HDMI Digital Output
- ☐ Red Laser Pointer
- ☐ Modular Design
- ☐ ESD Protected
- ☐ Stable Stands
- ☐ USB3.0 Converter
- ☐ Wide range of accessories
- ☐ Advanced Remote Controller
- ☐ Knob-Drive Remote Controller




Modular design of F30s allows easy configuration of your optical inspection system with a wide range of stands and mounts according to your needs.



F30s with Articulated Boom Stand



Double-arm Boom Stand

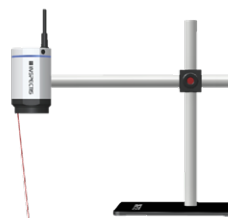

Single-arm Boom Stand



F30 with integrated stand

"Designed according to industry standards"

System Specifications

System:

Video Output	Full HD 1080/59.94p switchable to Full HD 1080/50p
Imager	3.2M pixels. 1920(H) x 1080(V) 16:9 format progressive scan CMOS image sensor
Pixel Size (resolution)	2.5 µm x 2.5 µm
I/O Ports	HDMI video output, 5P MiniDin Power input, 8P MiniDin External Remote Control
Onboard Control	Zoom, Iris, Auto Focus on/off, White Balance
External Remote Control (optional)	Zoom, Iris, Manual/Auto Brightness, Manual/Auto Focus, White Balance, Picture Quality Display Text, Crosshairs overlay, Image Rotate, Image Freeze, Memory Save/Recall (2 pcs), Foot-switch, Video Output Frame Rate
Illumination F30	Built-in LED matrix light. Optional LED RingLight.
Illumination F30s (optional)	LED RingLight, Single-Spot or Double-Spot LED matrix light with gooseneck
Storage Environment	-20° – +60 °C, 20 - 95% RH, non-condensing
Operating Environment	-5 ° – +40 °C, 20 - 95% RH, non-condensing
Size, Weight F30	395 x 240 x 360 mm including metal base, 4.5 kg. including metal base
Size, Weight F30s	φ 86.5 x 137mm exclusive cable. 1.1 kg. excluding stand

Optics:

Lens system	x30 Optical Zoom with Auto Focus and Iris Control
Working distance	228 mm (from lens to object at maximum zoom factor)
Magnification x1-x30 zoom factor	~ 1.9x - ~ 56x on 24" monitor (~ 170x with +10 macro lens)
Field of View at x1-x30 zoom factor	~ 270 x 151.9 mm to ~ 9.3 x 5.2 mm (~ 3.1 x 1.7mm with +10 lens)



Ordering numbers Inspection Systems:

Item name	Order No.		
Inspectis F30, Full HD 1080p Digital Macroscope, Integrated Stand, LED, HDMI	HD-012	Macro lens +5 diopter, 67mm (Achromatic)	HD-039
Inspectis F30s, Full HD 1080p Digital Macroscope, Modular, HDMI	HD-013	Macro lens +10 diopter, 67mm (Achromatic)	HD-123
With Antistatic housing	* -E	RingLight White LED with detachable diffuser	HD-024-B
With Laser Pointer	* -L	Polarization filter for Ring Light	HD-024-P-U
Upgrade with Hermetic Sealed Housing	* -S	RingLight UV 365-375nm wavelength LED	HD-025
Upgrade with 3m Extended Cable	* -EC	Single-Spot LED matrix light with gooseneck	HD-125
		Double-Spot LED matrix light with gooseneck	HD-126
		XY Floating Board, Antistatic	HD-028
		Full HD to USB3.0 Capture with Basics software	HD-031
		Full HD External Image Capture & Video Recorder	HD-140-F
		Single-arm Boom stand	HD-060
		Double-arm Boom Stand, Heavy Duty	HD-070
		Double-arm Boom Stand, Light	HD-071
		Articulated-arm Boom Stand	HD-080
		Track stand, Large Base with focusing, WD 35-410mm	HD-210-LB-E
		XY-travel Stand, In-Line	HD-170
		XY-travel Stand, Table-top	HD-180
		Rack-&-Pinion Fine Focusing module	HD-090
		Aluminium carrying case	HD-040

Ordering numbers optional Accessories:

Remote Camera Control with Advanced Adjustments for Inspectis Microscopes	HD-020
Foot switch for Remote Camera Control	HD-027
Remote Camera Control with Knob Drive	HD-200
Remote Camera Control with Knob Drive, Programable with 7 memories	HD-200-MEM
Macro lens +2 diopter, 67mm	HD-037
Macro lens +3 diopter, 67mm (Achromatic)	HD-038

Reseller: _____



ATEcare Service GmbH & Co. GmbH
Neufeldstr. 14
D-85232 Bergkirchen
+49 8131 318 575-120
info@atecare.com

Kontakt Schweiz:
ATEcare Alexander Hoertner
Friedhofweg 5
CH-9434 Au SG
+41 71 740 10 90
info@atecare.net

Internet:
www.atecare.com



Quality of Inspection

Inspectis Optical Systems
Stockholm, Sweden

www.inspect-is.com
info@inspect-is.com

